

HeelZup™ *Supreme*

Maximum-Lift Therapeutic Heel Elevating Cushion

One of the leading challenges for the wound care community is heel pressure ulceration. Heels are naturally susceptible to pressure ulcers and even the most advanced support surfaces prove to be ineffective at prevention and treatment. Existing modalities of heel pressure ulcer prevention and treatment have many drawbacks. Heel “boots” require frequent removal to inspect the skin integrity of the lower leg. Pillows, which slide out of place and compress very easily are also a high maintenance and unreliable solution to the problem. The HeelZup™ *Supreme* Therapeutic Heel Elevating Cushion provides the most practical solution to this complex problem:

- Extra thick construction provides higher elevation of the lower legs and heels.
- The 1.5” thick Visco-Elastic foam calf bed provides enhanced comfort and pressure reduction to the lower leg, from the achilles area to the knee.
- Raised side bolsters prevent legs from falling off of the sides of the cushion.
- Stays neatly in place and does not require constant maintenance and adjustment.
- 4-Way stretch cover allows the legs to immerse more deeply into the Visco-Elastic foam calf bed, while the No-Slip bottom material helps the cushion to stay neatly in place. No need for constant maintenance and adjustment.
- Manufactured of high-density, anti-microbial, hypo-allergenic and latex-free foam.



Phone: 866.736.8484
Fax: 941.359.6939
www.heelzup.com

Intensive Therapeutics, Inc.




U.S. Patent # 6,634,045

**P.O. Box 488
Sarasota, FL
34230**



The HeelZup® Supreme cushion provides higher elevation with enhanced comfort and pressure reduction, thanks to the Visco-Elastic Foam calf bed, (pictured below).

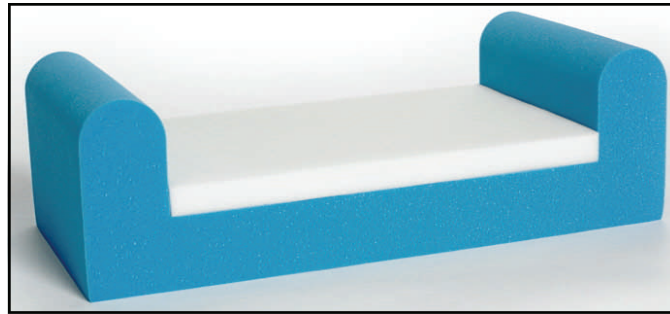

In this photo, the cover is removed to show the Visco-Elastic Foam calf bed. All HeelZup® Supreme cushions are reusable and fitted with a 4-Way Stretch Polyurethane cover.

Specifications:

HeelZup™ Supreme Reusable:

Manufactured of high-density, anti-microbial, hypo-allergenic foam. The standard width is 30" with 4" elevated leg bolsters on both sides. The calf bed will allow free lateral movement while the bolsters prevent the legs from falling off of the cushion. The foam base of the cushion is 3.5" thick, with a 1.5" layer of high-density, Visco-Elastic foam laminated onto the top surface. Total thickness is 5". A moisture proof, reusable 4-Way Stretch Polyurethane fabric cover with No-Slip bottom is sewn onto the cushion. The cover can be cleaned with standard cleaning procedures, using properly diluted solutions. HeelZup™ Supreme meets or exceeds the specifications in California Upholstery Provisions as stated in Technical Bulletin #117.

Ordering Information:

| Order Number | Description | Patient Height | Patient Weight |
|---------------------------------|--|----------------|---------------------|
| HZVE13RPS | HeelZup™ Supreme Reusable PetiteSoft™, 13" Depth | Up to 5'5" | 150 Pounds or Under |
| HZVE14R STANDARD SIZE | HeelZup™ Supreme Reusable, 14" Depth | Up to 6' Tall | Up to 300 Pounds |
| HZVE15R | HeelZup™ Supreme Reusable, 15" Depth | Over 6' Tall | Up to 300 Pounds |
| HZVE14RB | HeelZup™ Supreme Reusable Bariatric, 14" Depth | Up to 6' Tall | Over 300 Pounds |

Warranty:

This warranty applies only to defects in structure, workmanship and materials and only to damages arising from normal hospital usage on a non-prorated basis for a period of three (3) years for reusable models from the date of purchase. Damages arising from abnormal use such as those caused by needle punctures, burns, chemicals, negligent use or improper care or cleaning are excluded from coverage. Manufacturer will pay shipping charges incurred in connection with this warranty during the first year, thereafter, those charges are the facilities' or dealer's sole responsibility.